

WORLD INTELLECTUAL PROPERTY ORGANIZA International Burean



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(51) International Patent Classification 6:

G01N 33/50

(11) International Publication Number: WO 00/02045

A2

(43) International Publication Date: 13 January 2000 (13.01.00)

FP

(21) International Application Number: PCT/BE99/00087

(22) International Filing Date: 6 July 1999 (06.07.99)

22) International Fitting Date: 0 July 1999 (05,07.99)

6 July 1998 (06.07.98)

(71) Applicant (for all designated States except US): EURO-SCREEN S.A. [BE/BE]; Avenue des Becassines 7, B-1160 Brussels (BE).

(72) Inventors; and

(30) Priority Data:

98870149.6


(75) Inventors/Applicants (for US only): DUPRIEZ, Vincent [BE/BE]; Avenue de Roodebeek 229, B-1030 Brussels (BE). DETHEUX, Michel [BE/BE]; Rue Charles de Loupoigne 13/102, B-1348 Louvain-La-Neuve (BE). PARMENTIER, Marc [BE/BE]; Chaussée d'Uccle 304, B-1630 Brussels (BE).

(74) Agents: VAN MALDEREN, Eric et al.; Office Van Malderen, Place Reine Fabiola 6/1, B-1083 Brussels (BE). (81) Designated States: CA, IP, US, European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE).

Published

Without international search report and to be republished upon receipt of that report,

(54) Title: HIGH-THROUGHPUT SCREENING DIAGNOSTIC AND/OR DOSAGE METHOD OF AN AGONIST AND/OR AN ANTAGONIST FOR A CALCIUM-COUPLED RECEPTOR

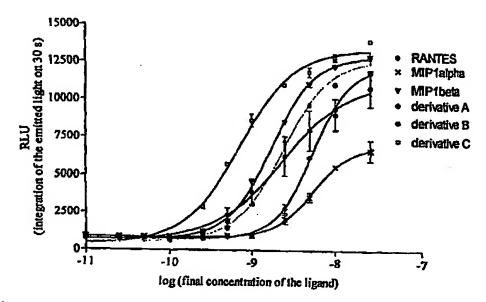

(57) Abstract

The present invention is related to a diagnostic and/or dosage method of an agonist and/or an antagonist and/or a modulator for a calcium-coupled receptor and/or channel and/or any other calcium-coupled protein, comprising the following successive steps: disposing the agonist and/or the antagonist and/or the modulator upon a solid support, incubating one or more cell(s) expressing aposequorin and said calcium-coupled receptor with coelenterazine in order to reconstitute an active acquorin in said cell(s); adding to said solid support one or more of said cells; and obtaining the measurement of an emitted light by said cell(s). The present invention is also related to the diagnostic and/or dosage device intended for the method according to the invention.

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
MA	Armenia	Pí	Finland	LT	Lithumia	SK	Slovakia
AT	Austria	FR	Prance	LÜ	Luxenbourg	SN	Scoegal
AU	Australia	GA	Gabon	LV	Latvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	Œ	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav	TM	Turkmenistan
BF	Burkina Faso	GR	Greece		Republic of Macedonia	TR	Turkey
BG	Bulgaria	HU	Hungary	ML	Mali	TT	Trinided and Tobago
BJ	Benin	IE	ireland	MN	Mongolia	UA.	Uknine
BR	Brazil	IL.	Israel	MR	Mauritania	UG	Uganda
BY	Belans	LS	Iceland	MW	Malgwi	US	United States of America
CA	Canada	IT	Italy	MX	Mexico	UZ	Uzbekistan
CF	Central African Republic	JР	Japan	NE	Niger	VN	Viet Nam
CG	Congo	KB	Keoya	NL	Netherlands	YU	Yogoslavia
CH	Switzerland	KG	Kyrgyzstan	NO	Norway	zw	Zimbabwe
a	Côte d'Ivoire	KP	Democratic People's	NZ	New Zealand		
СМ	Cameroon		Republic of Korea	PL	Poland		
CN	China	KR	Republic of Korea	PT	Portugal		
Cυ	Cuba	KZ	Kazakstan	RO	Romania		
cz	Czech Republic	LC	Saint Lucia	RU	Russian Federation		
DE	Germany	u	Liechtenstein	SD	Sudan		
DK	Denmark	ŁK	Sri Lanka	SE	Sweden		
ER	Estonia	LR	Liberia	SC	Singapore		

5

10

HIGH-THROUGHPUT SCREENING DIAGNOSTIC AND/OR DOSAGE METHOD OF AN AGONIST AND/OR AN ANTAGONIST FOR A CALCIUM-COUPLED RECEPTOR

Field of the invention

The present invention is related to a high15 throughput screening diagnostic and/or dosage method and
device of an agonist and/or an antagonist for a calciumcoupled receptor and the agonist and/or antagonist of said
calcium-coupled receptor identified by said method and
device.

20

Background of the invention and state of the art

A lot of G-protein-coupled receptors (GPCR) trigger, upon binding of an agonist, a transient increase in intracellular calcium concentration. This variation acts as an internal secondary messenger and is an important modulator of many physiological mechanisms (reviewed by Rink (1990), Tsunoda (1993) and by Santella & Carafoli (1997)). Measurement of intracellular calcium concentration in cells expressing a GPCR can thus be used to monitor the efficacy of activation of a GPCR by various compounds known - or suspected - to be a ligand for this GPCR.

The activation of other receptors such as ions channels may also induce an intracellular calcium concentration.

Changes in calcium concentration can be detected by several means and methods, like the use of fluorescent dyes (for example: fura-2, fluo-3, fluo-4 and indo-1).

However, Ca++ sensitive dyes have limitations. Activation of the dyes with an excitation beam requires complicated and expensive instruments and limits the use of the plastic labware such as microtiter plates.

Another method for intracellular calcium 10 concentration measurement is the use of cell lines overexpressing a GPCR and apoaequorin, such as described by Sheu et al. (1993). In this system, cells expressing apoaequorin are incubated with coelenterazine, which is the co-factor of aequorin. During this incubation, 15 coelenterazine enters the cell and conjugates with apoaequorin to form aequorin, which is the active form of the enzyme. Upon incubation of the cells with an agonist of the GPCR, intracellular calcium concentration increases. This increase leads to the activation of the catalytic 20 activity of aequorin, which oxidises coelenterazine and yields apoaequorin, coelenteramide, CO2 and light. Once the photon has been emitted, the complex must dissociate and apoaequorin must recombine with a new coelenterazine molecule to be able to emit light again. Thus, in this 25 system, measurement of light emission following agonist addition reflects its ability to activate the GPCR and thus to increase intracellular calcium concentration. Because light is emitted only during 20 to 30 seconds after activation of the GPCR, recording of the emitted light must 30 be performed during the few seconds following agonist addition to the cells. This flash-type signal is due to the fact that (1) intracellular calcium increase triggered by GPCR is only transient and (2) as mentioned earlier, after 10

oxidation of coelenterazine, apoaequorin must recombine with coelenterazine to be able to emit light again.

The Patent Application EP-0341477 teaches the expression of jellyfish photoprotein aequorin in a mammalian cell system by cloning gene pAQ440 specifying the biosynthesis of the aequorin into an expression vector plasmid of a mammalian cell system, subjecting the resulting plasmid to transfection and producing the photoprotein aequorin in the mammalian cell.

The Patent US-5,422,266 describes a gene encoding apoaequorin protein included in a vector capable of expressing the apoaequorin in E. coli.

The Patent US-5,714,666 describes mammalian cell lines or transgenic animals expressing apoaequorin and a receptor involved in the modulation of intracellular calcium. This document also describes a method of measuring intracellular calcium comprising adding coelenterazine cofactors to said mammalian cells expressing apoaequorin and measuring photoemission where emission of photons is indicative of intracellular calcium concentration.

However, the methods of the state of the art require firstly the spreading of cells from a mammalian cell line expressing apoaequorin on a solid support (for example a 96-well plate), secondly the addition of the coelenterazine cofactor upon the cells and incubation to reconstitute a functional aequorin, thirdly the preparation of the agent affecting a receptor involved in the modulation of intracellular calcium concentration, and its addition to the prepared cells, and finally the measurement of the photoemission.

Furthermore, as mentioned above, light is emitted only during 20 to 30 seconds after activation of the GPCR. Therefore, the recording of the emitted light



must be performed during the few seconds following agonist addition to the cells.

Therefore, the methods used in the state of the art are not adequate for a detection based upon highthroughput screening level, which usually need luminometer(s) and require the use of microtiter plates for the testing of thousands of compounds.

Aims of the invention

The present invention aims to provide a method and means, which do not present the drawbacks of the state of the art, for detecting biologically active substances, especially agonists and/or antagonists for calcium-coupled receptors.

A main aim of the present invention is to provide such method and means which allow the detection of biological active substances at a high-throughput scale, which could be adapted to specific recipients such as microtiter plates without requiring the modification of the high-throughput screening device.

Another aim of the present invention is to provide an easy and non-expensive method that could be easily automated.

25 Summary of the invention

The present invention is related to a high-throughput screening diagnostic and/or dosage method of an agonist and/or an antagonist for a "calcium-coupled" receptor, comprising the following successive steps:

- 30 disposing an agonist and/or an antagonist (preferably of a molecule) of said receptor upon a solid support,
 - incubating one or more cell(s) expressing apoaequorin and said "calcium-coupled" receptor with coelenterazine

in order to reconstitute an active aequorin by said cell(s).

- adding to said solid support one or more of said cells,
 and
- 5 obtaining the measurement of an emitted light by said cell(s).

The term "calcium-coupled" receptor means any receptor (such as a G-coupled receptor or an ion channel) whose activation (by an ion, a known or unknown agonist or antagonist molecule) may increase an intracellular calcium concentration in the cell comprising said receptor, preferably in its cytoplasmic membrane.

The terms "disposing... upon a solid support"
means the step of putting said compound into a recipient
such as a microtiter plate or any other solid support
without requiring any (covalent or other) fixing of said
compound to said solid support.

Advantageously, the solid support is a microtiter plate, preferably a 96-well microtiter plate.

Advantageously, the cell expressing appoacequorin and the calcium-coupled receptor is a cell expressing a G-coupled receptor and possibly one or more protein(s) intended to ensure a coupling of the receptor to the calcium pathway.

Preferably, said protein is selected from the group consisting of a natural Gal6 or Gal5 protein, a chimeric G-protein resulting from a fusion between two different G-proteins or a phospholipase Cβ2 protein.

The measurement of the emitted light is advantageously obtained with one or several luminometer(s), possibly equipped with several dispensers and measurement heads.

25



The present invention is also related to a high-throughput screening diagnostic and/or dosage device intended for the diagnostic and/or dosage method according to the invention, said device comprising the following 5 elements:

- a recipient, preferably a microtiter plate, more preferably a 96-well microtiter plate,
- a medium containing cell(s) expressing apoaequorin and a calcium-coupled receptor,
- 10 a medium containing coelenterazine, and
 - means (such as one or more luminometer(s) equipped with one or several dispensers and measurement heads) for detecting and possibly quantifying an emitted light by said cell(s).
- 15 Advantageously, the device according to the invention comprises means for the automatic performance of the successive steps of the diagnostic and/or dosage method according to the invention.

A last aspect of the invention is related to 20 the agonist and/or the antagonist of a calcium-coupled receptor identified by the method or the device according to the invention.

The present invention will be described in details in the following non-limiting examples, in 25 reference to the enclosed figures.

Brief description of the drawings

Figure 1 shows a series curves representing the intensity of the emitted light by cells as a function of time for each well of a 96-well plate injected with cells expressing the CCR5 receptor, apoaequorin and Ga16. The scaling is the same for all the graphs. Recording of

5

the signal was performed for 30 seconds. Ligand concentrations are increasing from column 1 towards column 12. All measurements were performed in duplicate: lines A and B: ligand is RANTES; lines C and D: ligand is MIP-1 α ; lines E and F: ligand is MIP-1 β ; lines G and H: ligand is derivative A of RANTES.

- Figure 2 represents the dose-response curve for different agonists of the CCR5 receptor which represent RLU (integration of emitted light on 30 seconds) according to the logarithm of the final concentration of the ligand.
- Figure 3 represents the dose-response curve for different agonists of the 5HT-2B receptor.
 - Figure 4 represents the dose-response curve of the light emitted from K562 cells expressing CCR3 and aequorin as a response to the activation of the receptor by eotaxin.
- 20 Figure 5 represents the dose-response curves for cells expressing aequorin and Gα16 and (panel A) the orexin 1 receptor or (panel B) the orexin 2 receptor.
- Figure 6 represents the dose-response curve for different antagonists of the 5HT-2B receptor.

Description of a preferred embodiment of the present invention

Detection of agonistic activities by means of
mammalian cell lines expressing apoaequorin and a GPCR
requires the measurement of the emitted light to be
performed just after placing the cells in contact with the
potential agonist. This can easily be measured at low

throughput using a single-tube luminometer. However, up to now, this biological system could not be used at a high-throughput scale. Indeed:

- the necessity to measure light just after placing the 5 cells in contact with the agonist to be tested compels to use a luminometer equipped with a build-in dispenser. For example, due to the short duration of light emission, it is impossible to inject the drugs to be tested on the cells placed in the 96 wells while 10 plate is outside the luminometer and subsequently record emitted light with the plate into the luminometer. Even if the plate could rapidly (i.e. less than 15 seconds) be placed luminometer after injection of the drugs to be tested, 15 current apparatus do not allow the measurement of light from the 96 wells before the extinction of the flash signal of aequorin, as these luminometers are not equipped with 96 detectors.
- 20 allow to inject a single solution into the 96 wells, making it impossible to inject a different drug in each well. Moreover, the washing of the dispenser before each measurement, for the injection of another drug in the next well, is time-consuming and thus is not suitable for the high-throughput scale. The same problem occurs with devices equipped with 6 dispensers (e.g. the "Microbeta Jet" from EG&G Wallac) as the 6 dispensers only deliver a single solution.

The present invention provides a method for 30 performing high-throughput screening of drugs binding to GPCR by the use of mammalian cell lines expressing apoaequorin and a GPCR and by the use of a conventional luminometer. Following this method, the solutions to be tested for (ant)agonistic activities are placed in the

10

wells of a 96-well plate. Cells expressing apoaequorin and a GPCR are detached from the culture plate (or collected cultures) and are incubated suspension coelenterazine to reconstitute active aequorin. These are 5 then maintained in suspension with a magnetic stirrer and the cell suspension is injected, well by well, on the solutions of potential agonist to be tested. Light emission is then recorded for 1 (alternatively up to 30 or more) seconds. This method, by injecting the same cell suspension in each of the 96 wells, avoids the need of washing the dispenser between each measurement and allows to perform 96 measurements of agonist-induced aequorin light emission in 15 minutes or less with a single dispenser luminometer. Alternatively, it allows to perform 96 measurements of 15 agonist-induced aequorin light emission in 2 minutes or less with a luminometer equipped with 6 dispensers and measurement heads (e.g. with the "Microbeta Jet" from EG&G Wallac).

This method thus allows to perform highthroughput screening (10 000 samples/day) with mammalian cell lines expressing apoaequorin and a GPCR and by the use of a conventional luminometer. This reduces the screening time and the amount of drugs needed for each measurement.

perform a also allows to system 25 functional screening with very few (down to 5000 or less) cells per measurement.

The injection of the cells into the wells containing the agonists did not increase the background of the measurement (that could for example have originated from cell breakage, releasing aequorin molecules from the into the culture medium, where the concentration would have triggered the emission of light from aequorin). A signal-to-noise ratio above 50 was commonly obtained with this system of cell injection.



The method according to the invention is suitable for performing high-throughput analysis of GPCR or other calcium-coupled-receptor stimulation by known or agonists by means of cells expressing the 5 receptor and apoaequorin. These cells apoaequorin in the cytoplasm, as described by Sheu et al. (1993) or Button and Brownstein (1993) or may express apoaequorin in the mitochondria, by means of the addition of a mitochondrial targeting sequence to the aequorin, as 10 used by Stables et al. (1997) or in any other part of the cell. These cells may also express proteins intended to ensure coupling of the over-expressed receptor to the calcium pathway. These may be the natural Gal6 or Gal5 proteins (Milligan et al., 1996), chimeric G proteins 15 resulting from a fusion between two different G proteins (Komatsuzaki et al., 1997), phospholipase C-β2 (Park et al., 1992), or any other "universal coupling" protein. Once the cells have been prepared and loaded with coelenterazine, they can be used for several hours (at **20** least 9 hours). The load in coelenterazine and the intensity of the light emitted by the cells upon agonist stimulation lasts is stable for this period of time.

Examples

25 Example 1

A CHO cell line expressing the chemokine CCR5 receptor, the Gα16 coupling protein and apoaequorin was
established. Cells were cultivated as a monolayer in
HAM'sF12 medium containing 10% Foetal bovine serum (FBS).

30 On the day of the experiment, the culture medium was
removed and cells were incubated for 5 min at room
temperature in PBS-EDTA (phosphate buffered saline solution
without calcium, supplemented with 5 mM EDTA). Cells were

detached from the culture vessel by shaking the culture plate by hand and by pipetting up and down. Cells were centrifuged and the supernatant was removed to eliminate the EDTA; the pellet was resuspended in HAM's F12 culture medium without FBS and with 0.1% Bovine Serum Albumin. Cells were counted by means of a Thomas cell, were centrifuged again and were resuspended in HAM's F12 culture medium without FBS and with 0.1% Bovine Serum Albumin at a concentration of 5.106 cells/ml. Coelenterazine (or a derivative of it, e.g. Coelenterazines f, h, n, cp or hcp, from Molecular Probes Inc.) at 500 µM in methanol was added to the cell suspension at a final concentration of 5 µM. The cell suspension was then stored in the dark at room temperature for 3 to 5 h, with shaking every 15 to 30 min to maintain the cells in suspension.

Series of dilutions of known ligands were prepared in HAM's F12 culture medium without FBS and with 0.1% Bovine Serum Albumin and 50 µl of each of these solutions were placed in the wells of a 96-well plate. The cell suspension was diluted 5 times with medium HAM'sF12 without FBS and with 0.1% Bovine Serum Albumin and was placed in a glass or plastic container protected from light by wrapping it with aluminium foil. A magnetic stirring bar was added to the suspension and a magnetic stirrer was used at low speed (1 to 5 rounds per second) to maintain the cells in an homogenous suspension. The magnetic stirring bar was equipped with a ring to avoid crushing the cells, and the subsequent release of aequorin in the culture medium. Alternatively, a culture vessel equipped for culture of cells in suspension may be used.

One uses the EG&G Wallac's MicroLumat-Plus microplate luminometer, which allows injection and direct subsequent recording of the emitted light from each well of



a 96-well plate. The end of the entrance tube of the dispenser was placed at the bottom of the cell suspension and the dispenser was washed with 3 times the dead volume of the apparatus so that all the volume of the tube and 5 pumps was filled with the cell suspension. The 96-well plate containing the solutions of agonists was inserted into the luminometer. Then, for each well, 50 μ l of the cell suspension (i.e. 100 000 cells) was dispensed into the well (at the lowest injection speed (0.4 s) to 10 prevent cell breakage that would release aequorin into the culture medium) and the emitted light was immediately recorded during 30 seconds. After reading the first well, cells were injected into the next well and emitted light was recorded, etc. For each plate, a series of curves 15 representing the intensity of the emitted light as a function of time for each well was displayed (figure 1). The intensity of the emitted light was integrated over 30 s using the Winglow software provided with the luminometer, yielding for each well one value representative of the 20 emitted light and hence of the stimulation of the CCR-5 receptor by the agonist present in the well. These values can be plotted against the logarithm of the ligand concentration to generate dose-response curve as shown in figure 2. These allow the determination of half-maximal 25 response doses (EC50) for each ligand. For the generation of these data, 288 measurements were performed in less than 3 hours.

Example 2

A CHO cell line expressing the serotonin 5HT-2B receptor, the Gal6 coupling protein and apoaequorin was established. Cells were treated as described in Example 1 and were dispensed after dilution (100 μ l/well, WO 00/02045

corresponding to 50 000 cells) on 100 μ l of solutions of known agonists for this receptor. The emitted light was recorded during 20 s for each well. Dose-response curves obtained for different agonists are shown in figure 3. For the generation of these data, 144 measurements were performed in less than 1 hour.

13

Example 3

K562 cells expressing the chemokine CCR3

10 receptor were transfected by a plasmid for the expression of aequorin and the Gα16 coupling protein. Cells stably transfected were selected for 2 weeks with the antibiotic Zeocin. These cells were cultured in suspension in DMEM culture medium containing 10% FBS. They were centrifuged and the pellet was used as described in example 1 to perform aequorin measurements. A dose-response curve with eotaxin and MCP-4 generated by this method is described in figure 4.

20 Example 4

A CHO cell line expressing the orexin 1 or orexin 2 receptor, the Gα16 coupling protein and apoaequorin was established. Cells were treated as described in Example 1 and were dispensed after dilution (100 μl/well, corresponding to 25 000 cells) on 50 μl of solutions of known agonists for these receptor. The emitted light was recorded during 20 s for each well. A doseresponse curve obtained with this method is showed in figure 5 A.

30

<u>Example 5</u>

 $$\tt A$$ CHO cell line expressing the canabinoid CB1 receptor, the G\alpha16 coupling protein and apoaequorin was



established. Cells were treated as described in Example 1 and were dispensed after dilution (100 μl/well, corresponding to 50 000 cells) on 100 μl of solutions of a known agonist for this receptor. The emitted light was recorded during 20 s for each well. A dose-response curve obtained with this method is showed in figure 5 B.

For the measure of antagonistic activities. cells were injected on series of dilutions of antagonists. At this stage, the emitted light was recorded to check that 10 the potential antagonists had no agonistic activity. The cells were incubated with the antagonists for 30 min. A solution of agonist (alpha-methyl-5HT) of the receptor was then injected on the cells and the emitted light was immediately recorded for each well. The emitted 15 light was plotted as a function of the logarithm of the concentration to yield the graph antagonist Increasing antagonist concentrations result in a decreasing light emission upon agonist addition. An agonist of another receptor expressed by the cells (usually ATP, acting at P2 20 receptors) can then be injected on the mixture antagonist, cells and agonists as a control to check that the cells still have active aequorin up to that moment of the experiment. For example, cytotoxic compounds that increase the intracellular calcium concentration can make 25 the aequorin consume all the coelenterazine present in the cell. Such a cytotoxic compound will be detected by the absence of signal upon ATP injection (Fig. 6).


REFERENCES

- Button, D. & Brownstein, M., Cell Calcium 14, pp. 663-671 (1993)
- Komatsuzaki, K. et al., FEBS Lett. 406, pp. 165-170
 5 (1997)
- Milligan, G. et al., TIPS 17, pp. 235-237 (1996)
 - Park, D. et al, J. Biol. Chem. 267, pp. 16048-16055 (1992).
 - Rink, T.J., FEBS Lett. 268, pp. 381-385 (1990)
- - Sheu, Y.-A. et al, Analytical Biochemistry 209, pp. 343-347 (1993)
- Stables, J. et al., Analytical Biochemistry 252, pp. 15 115-126 (1997).
 - Tsunoda, Y., Biochim. Biophys. Acta 1154, pp. 105-156 (1993)



CLAIMS

- Diagnostic and/or dosage method of an agonist and/or an antagonist or any modulator for a calcium-coupled receptor or a calcium-coupled channel or any other calcium-coupled protein, comprising the following successive steps:
 - disposing the agonist and/or the antagonist upon a solid support,
- incubating one or more cell(s) expressing apoaequorin or any other related protein and said calcium-coupled receptor with coelenterazine or any other cofactor of a calcium-sensitive protein in order to reconstitute an active aequorin by said cell(s),
- adding to said solid support one or more of said cells,
 and
 - obtaining the measurement of an emitted light by said cell(s).
 - 2. Method according to claim 1, wherein the solid support is a microtiter plate.
- 20 3. Method according to claim 2, characterised in that the microtiter plate is a 96-well microtiter plate, or a 384-well plate, or a 1536-well-plate or any other format.
- 4. Method according to any one of the 25 preceding claims, characterised in that the cell expresses apoaequorin in the cytoplasm or in the mitochondria or in any other part of the cell.
- 5. Method according to any one of the preceding claims, wherein the cell expressing a calcium-30 coupled receptor is a cell expressing an endogenous or recombinant G-protein-coupled receptor and/or a cell which expresses proteins intended to ensure a coupling of the

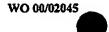



analysed receptor (endogenous or overexpressed) to the calcium pathway.


- 6. Method according to claim 5, wherein said protein is selected from the group consisting of natural 5 Gα16 or Gα15 protein, chimeric G-protein resulting from a fusion between two different G-proteins or phospholipase Cβ2 protein or any other coupling protein or chemical.
- 7. Method according to any one of the preceding claims, characterised in that the measurement of the emitted light is obtained with one or more luminometer(s), advantageously equipped with several dispensers and measurement heads.
- 8. High-throughput screening diagnostic and/or dosage device intended for the high-throughput screening diagnostic and/or dosage method according to any one of the preceding claims, comprising the following elements:
 - a microtiter plate, preferably a 96-well microtiter plate,
- 20 a medium containing cell(s) expressing apoaequorin and a calcium-coupled receptor,
 - a medium containing coelenterazine, and
 - means for detecting an emitted light by said cell(s).
- 9. Device according to claim 8, comprising 25 means for the automatic performance of the successive steps of the diagnostic and/or dosage method according to any one of the claims 1 to 7.
- 10. Agonist or antagonist of a receptor identified by the method according to any one of the claims 30 1 to 8.

1/6

FIG.1

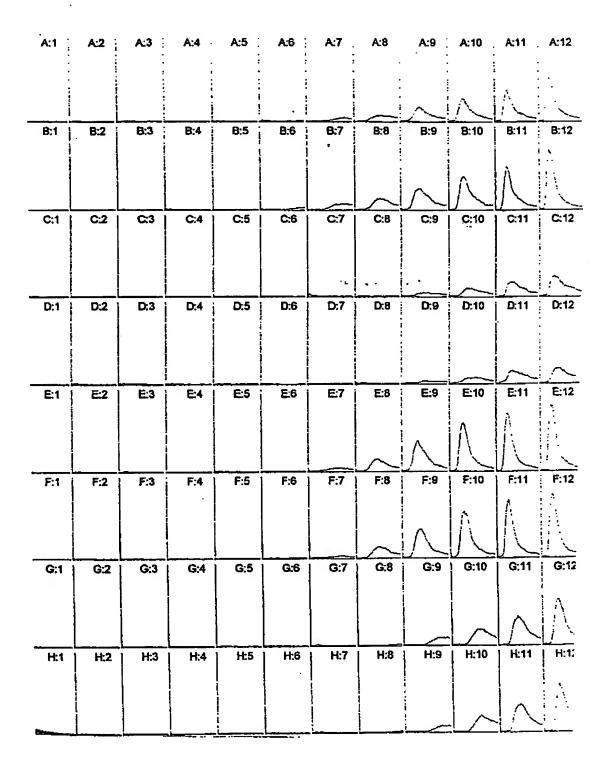

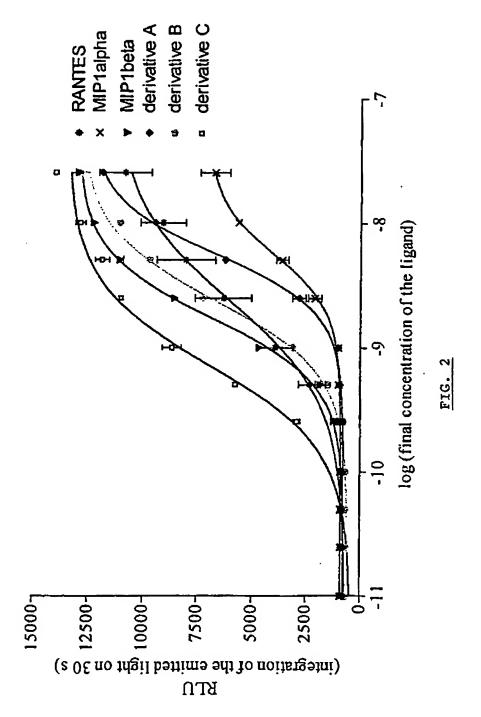

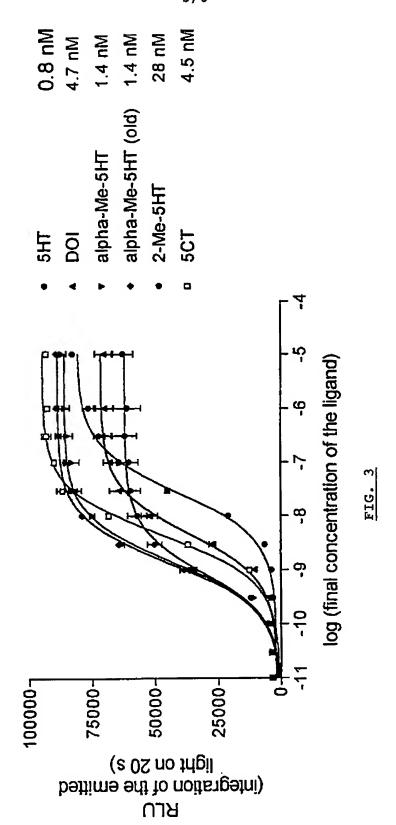

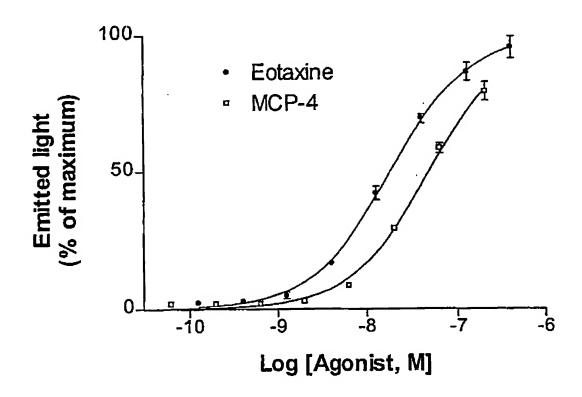


FIG. 4

WO 00/02045

A.

Dose response curve Mix Aequorin-Orexin 1 R

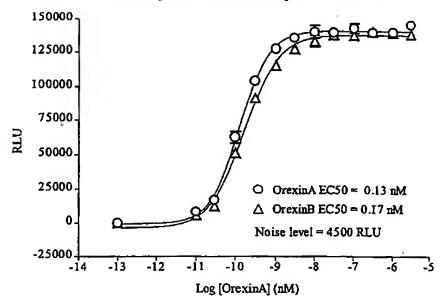


FIG. 5A

B.

Dose response curve Mix Aequorin-Orexin 2 R


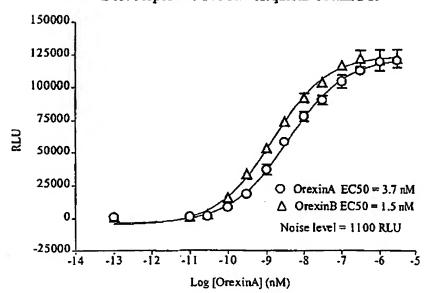


FIG. 5B

7.

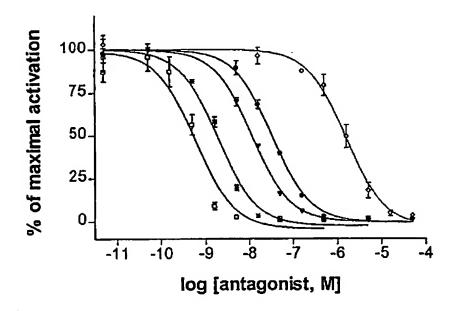

- Methysergide
- Methiothepin
- mesulergine
- Mianserin
- Ketanserine

FIG. 6

(19) World Intellectual Property Organization International Bureau



(43) International Publication Date 13 January 2000 (13.01.2000)

PCT

(10) International Publication Number WO 00/02045 A3

(51) International Patent Classification6: 33/566, 33/543

G01N 33/50.

Loupoigne 13/102, B-1348 Louvain-La-Neuve (BE). PARMENTIER, Mare [BE/BE]; Chaussée d'Uccle 304, B-1630 Brussels (BE).

(74) Agents: VAN MALDEREN, Eric et al.; Office Van Malderen, Place Reine Fabiola 6/1, B-1083 Brussels (BE).

(84) Designated States (regional): European patent (AT, BE,

CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC,

(81) Designated States (national): CA, JP, US.

- (21) International Application Number:
 - PCT/BE99/00087
- (22) International Filing Date:
- 6 July 1999 (06.07.1999)
- (25) Filing Language:

English

(26) Publication Language:

English

(30) Priority Data:

98870149.6

6 July 1998 (06.07.1998)

(71) Applicant (for all designated States except US): EURO-

Published:


NL, PT, SE).

- with international search report
- SCREEN S.A. [BE/BE]; Route de Lennik 802, B-1160 Brussels (BE).
- (88) Date of publication of the international search report: 8 November 2001

(72) Inventors; and (75) Inventors/Applicants (for US only): DUPRIEZ, Vincent [BE/BE]; Avenue de Roodebeek 229, B-1030 Brussels (BE). DETHEUX, Michel [BE/BE]; Rue Charles de

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

(54) Title: BIOLUMINESCENT ASSAY FOR AGONISTS OR ANTAGONISTS OF A CALCIUM-COUPLED RECEPTOR

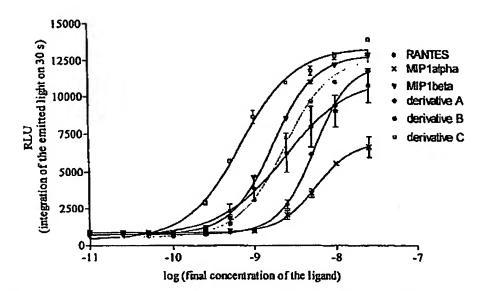

(57) Abstract: The present invention is related to a diagnostic and/or dosage method of an agonist and/or an antagonist and/or a modulator for a calcium-coupled receptor and/or channel and/or any other calcium-coupled protein, comprising the following successive steps: disposing the agonist and/or the antagonist and/or the modulator upon a solid support; incubating one or more cell(s) expressing apoaequorin and said calcium-coupled receptor with coelenterazine in order to reconstitute an active aequorin in said cell(s); adding to said solid support one or more of said cells; and obtaining the measurement of an emitted light by said cell(s). The present invention is also related to the diagnostic and/or dosage device intended for the method according to the invention.

INTERNATIONAL SEARCH REPORT

tr. attenut Application No

		PC1/BE 99	/0008/		
A CLASSI IPC 7	FICATION OF SUBJECT MATTER G01N33/50 G01N33/566 G01N3	3/543	-		
According to	o International Paterni Classification (IPC) or to both netional class	ssification and IPC			
B. FRELDS	SEARCHED				
Minimum de IPC 7	commentation searched (classification system followed by classification sy	fication symbots)			
	fion searched other than minimum documentation to the extent t				
Bectronic d	ata base consulted during the international search (name of dat	ia base and, where practical, search terms used)		
C. DOCUM	ENTS CONSIDERED TO BE RELEVANT				
Calegory *	Citation of document, with indication, where appropriate, of the	ne relevant passagės	Relevant to claim No.		
X	CHEMICAL ABSTRACTS, vol. 127, 18 December 1997 (1997-12-08) Columbus, Ohio, US; abstract no. 314897, XP002121177 abstract	no. 23,	1-10		
	& E. KNIGHT ET AL: "A biolumi: assay for agonist activity at any G-protein-coupled receptor ANALYTICAL BIOCHEMISTRY, vol. 252, no. 1, 1997, pages 1: New York NY USA				
Furt	her documents are listed in the continuation of box C.	Patent family members are listed	in annex.		
Special ca	siegories of clied documents :	and before the second subset of the second			
"A" docume consid	eril defining the general state of the art which is not tered to be of particular relevance document but published on or after the international	or priority date and not in conflict with cited to understand the principle or the invention	"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention "X" document of particular relevance; the claimed invention		
which citalio	ent which may throw doubts on priority claim(s) or is caled to establish the publication date of another n or other special reason (as specified)	cannot be considered novel or earmot involve an inventive step when the do "Y" document of particular relevance; the c cannot be considered to involve an in-	cannot be considered novel or carmot be considered to involve an inventive step when the document is taken alone "Y" document of particular relevance; the claimed invention carnot be considered to involve an inventive step when the		
other)	ent referring to an oral disclosure, use, exhibition or means ent published prior to the international filing date but han the priority date claimed	ments, such combination being obviou in the art.	document is combined with one or more other such docu- ments, such combination being obvious to a person skilled in the art. *&* document member of the same patent family		
	actual completion of the International search		Date of mailing of the International search report		
	November 1999	11/11/1999			
Name and r	mailing address of the ISA European Palent Office, P.B. 5818 Patentiaen 2 NL – 2280 HV Riswijk	Authorized officer	Authorized officer		
	Tel. (+31-70) 340-2040, Tx. 31 651 epo nl, Fex: (+31-70) 340-3016	Yan Bohemen, C	∀an Bohemen, C		

Form PCT/ISA/210 (second sheet) (July 1992)

1

This Page is Inserted by IFW Indexing and Scanning Operations and is not part of the Official Record

BEST AVAILABLE IMAGES

Defective images within this document are accurate representations of the original documents submitted by the applicant.

Defects in the images include but are not limited to the items checked:

□ BLACK BORDERS
 □ IMAGE CUT OFF AT TOP, BOTTOM OR SIDES
 □ FADED TEXT OR DRAWING
 □ BLURRED OR ILLEGIBLE TEXT OR DRAWING
 □ SKEWED/SLANTED IMAGES
 □ COLOR OR BLACK AND WHITE PHOTOGRAPHS
 □ GRAY SCALE DOCUMENTS
 □ LINES OR MARKS ON ORIGINAL DOCUMENT

IMAGES ARE BEST AVAILABLE COPY.

☐ OTHER:

As rescanning these documents will not correct the image problems checked, please do not report these problems to the IFW Image Problem Mailbox.

☐ REFERENCE(S) OR EXHIBIT(S) SUBMITTED ARE POOR QUALITY